

Instrument Cassettes


Save **time**. Improve **reliability**.

In contrast to trays, cassettes only need to be set up *once*; post-procedure, the cassette is processed as a complete, sealed unit. Eliminating the breakdown/setup cycle in this manner saves an enormous amount of time and effort associated with every procedure.

Cassettes also remove a critical margin of error, because instruments are never removed from the cassette outside the operator – *and an instrument is assigned to every slot*. This allows the practitioner to determine, at a glance, if something is missing.

High flow-through design and stainless steel construction with full-width hinges make these units both efficient and durable.

Thompson 5 Cassette (item M2G1)

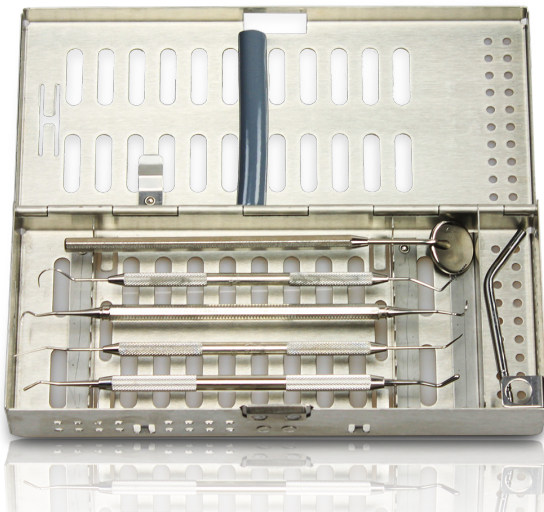

Exam, Ortho, Pedo-hygiene & Quick-Check Cassette

This setup covers preliminary exam and diagnostic instruments.


Configuration: 5 Instrument Slots,
plus Side Compartment

Dimensions: 8" x 3-3/16" x 1-1/8" (closed)

Typical Setup: Mirror, Explorer, Scaler, Perio Probe,
Composite Instrument, and A/W Syringe Tip

Also available in a 7-slot Slimline version (M2A1)

9 Instrument Cassette (item M2B2)

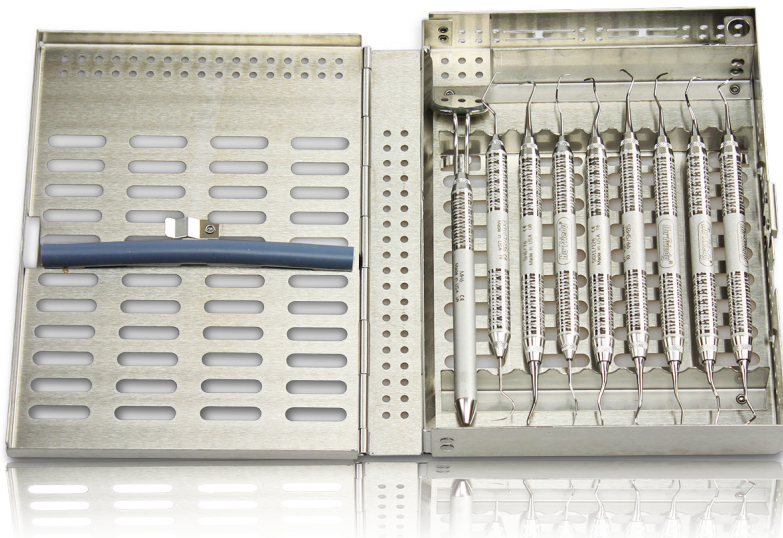

Hygiene Cassette

This cassette is the only unit you need for *any* hygiene procedure (Prophy, Perio, etc). Don't complicate matters by over-specializing – always think consolidation.

Configuration: 9 Instrument Slots
plus Side Compartment


Dimensions: 8" x 6" x 1-3/16" (closed)

Typical Setup: Mirror, Explorer, Perio Probe,
Posterior Explorer, Anterior Sickle,
Gracey 1/2, Universal 5/6, Gracey 11/12,
Gracey 13/14, and A/W Syringe Tip

*All Instruments are shown for display purposes only and are not included in any Cassette purchase.

Instrument Cassettes

14 Instrument Cassette (item M2C1)

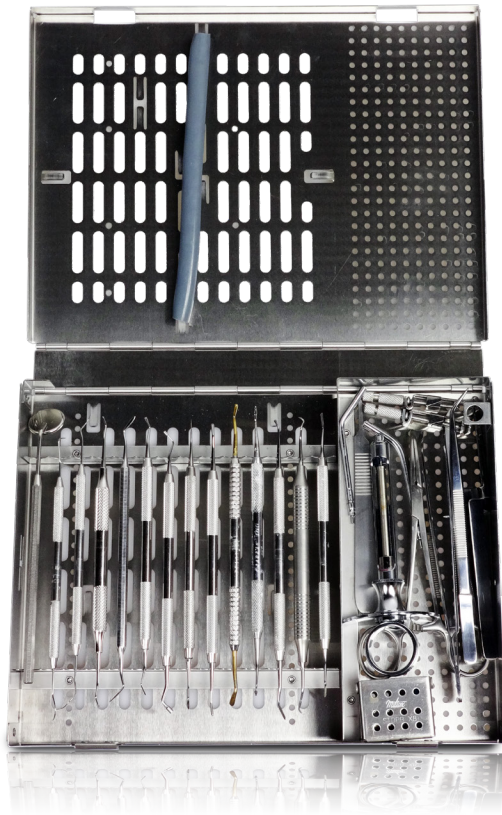

Restorative/Operative/Crown & Bridge Cassette

This unit is your bread and butter cassette and is used for all restorative dentistry.

Configuration: 14 Instrument Slots plus Side Compartment

Dimensions: 8" x 10-7/8" x 1-3/16" (closed)

Typical Restorative Setup:

Slots - Mirror, Perio Probe, Cowhorn Explorer, Spoon Excavator, Hollenback Carver, Cleoid/Discoïd Carver, Amalgam Condenser Plugger, Interproximal Carver, Small Plugger, Pear-shaped Burnisher, Finishing Burnisher, Composite Instrument, Universal Posterior Sickle Scaler, Packing Instrument

Side Compartment - Anesthetic Syringe, Articulating Paper Holder, 2 Matrix Band Holders, A/W Syringe Tip, Hemostat, 12-Hook Scalpel, 2 Cotton Forceps

Cubed Container - Air Abrasion Tip

Also available in a 16-slot version (M2C2)

Surgical Cassette (item M2E1)



Medium size shown above. Surgical Cassettes are also available in a larger size (M2F1), and smaller size (M2D1).

Ideal for large, bulky surgical instruments. Surgeries also tend to happen less frequently in most practices – making them great candidates for a “set it and forget it” cassette.

Configuration: 6 Instrument Slots plus 2 Compartments and Separator.

Dimensions: 10" x 5-1/2" x 4" (closed)

Typical Restorative Setup:

Slots - Mirror, Explorer, Periosteal, Luxating Spade Elevator, East/West Elevators

Left Compartment - Syringe, Hemostat, Scissors

Main Right Compartment - Curette, Small Straight Elevator, Large Straight Elevators, Bone File, Upper Forceps, Lower Forceps, Anterior Forceps, Closed Beak Forceps, Hemostat

Small Right Compartment - Tissue Forceps, A/W Syringe Tip, Surgical Scalpel, Surgical Scissors