

TM

ergonomic productsTM
Performance Without Compromise[®]



Lighting Solutions

The following pages excerpted from Ergonomic Products catalog: 950-075revM

[Click here to download full catalog](#)

www.ergonomic-products.com



All-in-one
functionality that
is **brilliant**
in every sense
of the word.



Aurora LS™

A truly exceptional lighting solution addresses both clinical needs and patient comfort – which is why Ergonomic Products has developed the TBLco brand (“The Best Light”), and bring you the Aurora LS.

The Aurora LS is the industry’s only fully integrated light incorporating full room, spot, and intraoral lighting with patient media. This unit condenses all illumination needs into a single, elegant fixture.

The Aurora’s dimmable room strips cast light both up and to the sides for an even, reflected coverage, making patients more comfortable and reducing clinician eye strain.

Using all LED lights, the Aurora LS is significantly more energy efficient than comparable products, generates less heat, and conforms to the most stringent state regulations.



Standard Features

- All LED lighting
- 8000-30000 lux Intraoral Light
- Built-in touchless On/Off control
- Tuneable wave-length allows intraoral use with a range of curable composites
- Continuous 360° Light-arm rotation
- Dimmer-controlled, reflected General Lighting for an improved patient experience
- 2 high-intensity Task Lights, with independent direction and focus adjustment
- Industry exclusive igus® cable management prevents cord binding and wear
- Adjustable VESA 100 Monitor Bracket (monitor not included)





Dimensions & Specifications

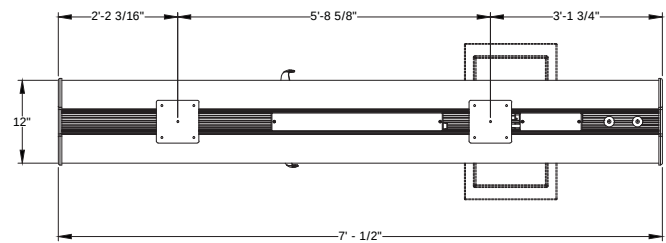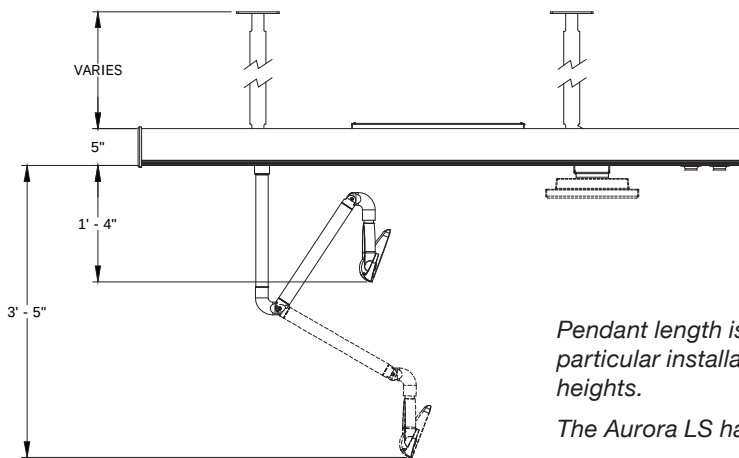

Pendant length is dependent on ceiling height and is fabricated to suit your particular installation. You will need to provide structural and finished ceiling heights.

The Aurora LS has a minimum ceiling height requirement of 9 feet.

⚡ The Aurora LS Track Light has 4 individual circuits as listed below:

LED General Room Lighting:

120v, 60hz
198w*
Color Temp: 5000k
Dimmable: 0-10vΩ

LED Oral Cavity Light:

120v, 60hz
10w
Color Temp: 3500k and 5500k
Dimmable: No

LED Task Lighting (2):


120v, 60hz
8w, 3-4vdc, 680mA
Color Temp: 4000k
Dimmable: No

TV Outlet:

120v, 60hz
300w


*reducible to 120w for CA 24 compliance

Ideal Room Location

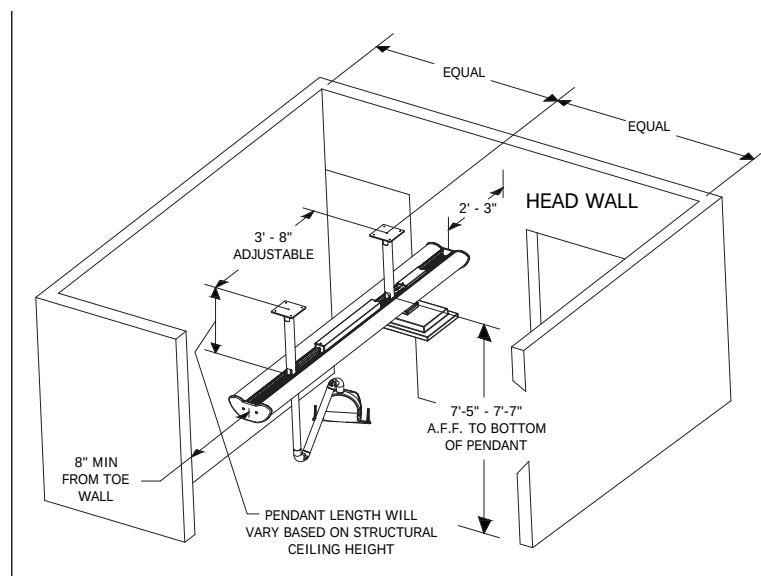

The Aurora LS is intended for center-line mounting in the operatory.

As a reflected-light system, it requires pendant mounting. Ideal height from floor to the bottom of the pendant plate is between 7'5" and 7'7".

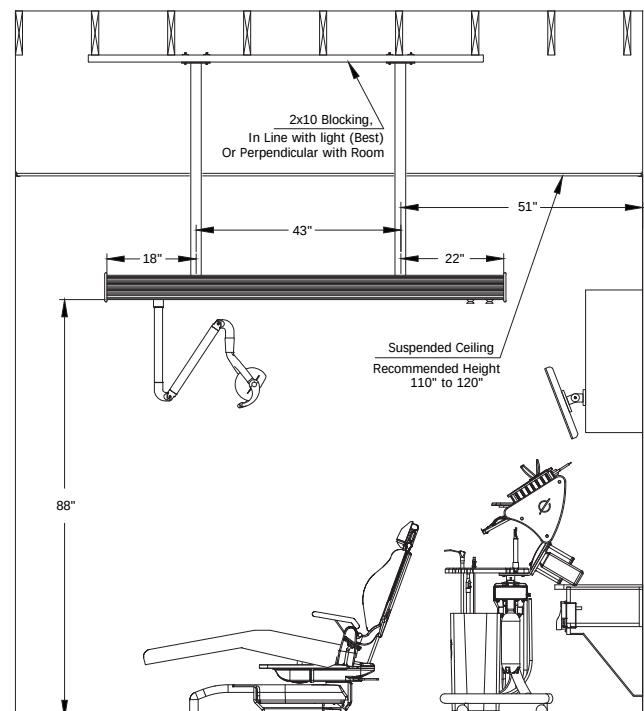

FastTrack LS™ (Item L2A1)



FastTrack LS™

The FastTrack LS is an extremely cost-effective track-light that provides high-quality features at a very attractive price.

Designed with industry-leading, tuneable wave-length LED output, the FastTrack LS provides enhanced visual acuity while reducing heat generation and lowering energy costs.

And our touchless light controls make asepsis management as simple as a wave of the hand!

With its milled aluminum trolley, single-piece extrusion body, 360° arm-rotation and an advanced igus® cable management system, this light is built to out-perform – and out-last – the competition.

Easy to install. Built to last.

Our intraoral light head, on a **simple, compact & reliable** track for **exceptional lighting** at an **affordable price**.



**Don't be left
in the dark!**

The **FastTrack LS** installs in half the time of other lighting systems, making it an exceptional replacement for older, outdated lights. Scan the QR code above to see how quickly a single installer can get you up and running!

Standard Features

- All LED lighting
- 8000-30000 lux Intraoral light
- Built-in touchless On/Off control
- Tune-able wave-length allows intraoral use with a range of curable composites
- 360° Light-arm rotation
- Single-piece aluminum extrusion Body
- Milled billet aluminum Trolley
- Industry exclusive igus® cable management prevents cord binding and wear



ergonomic products
Performance Without Compromise™



www.ergonomic-products.com • 866-374-6487

Excerpted from EP Catalog:

950-075revM



Dimensions & Specifications

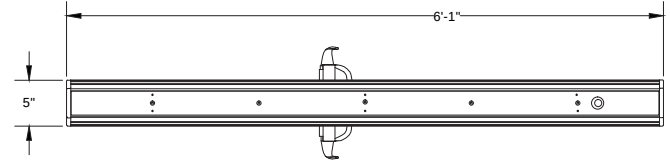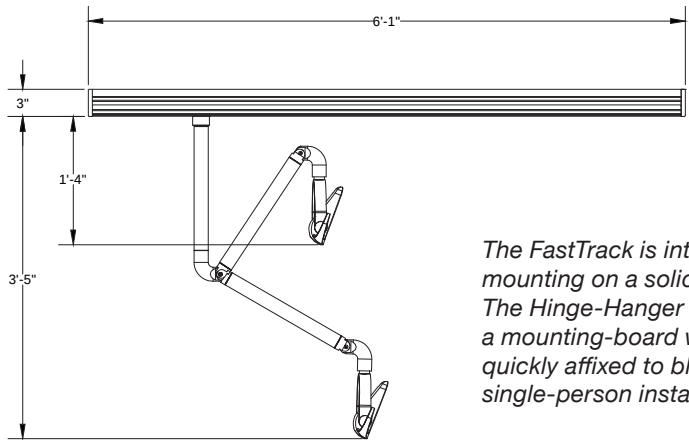


The FastTrack is intended for mounting on a solid, flat ceiling. The Hinge-Hanger system includes a mounting-board which may be quickly affixed to blocking allowing single-person installation.

Electrical Specifications

120v, 60hz
10w
Color Temp: 3500k & 5500k

Ideal Room Location

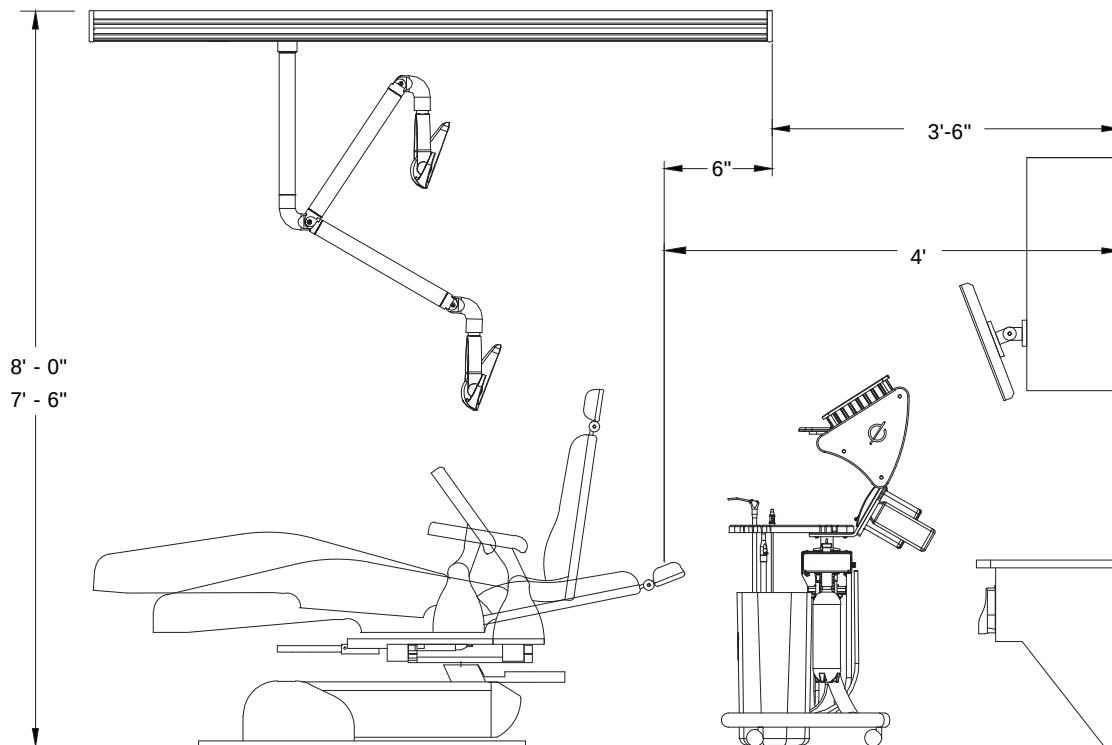

The FastTrack is intended for center-line mounting in the operator.

The layout above is represents a typical EP Workstation-equipped operator. Headwall-to-Track distances may vary in other environments. If so, position light based on reclined patient chair head-rest as indicated.

FastTrack LS™ PLUS (FastTrack upgrade path)



The same durability and **ease-of-installation** features offered on our FastTrack... but with **enhanced clinical function** and an **improved patient experience**.

TBL_{co}

FastTrack LS™ PLUS

Take your operatory lighting to the next level.

The FastTrack LS PLUS incorporates 2, high-intensity LED task-lights operated by a wireless remote.

Add an optional monitor mount, and this next-generation light really shines, while remaining simple to install, and *still* competitively priced.

Standard Features

- All LED lighting
- 8000-30000 lux Intraoral light
- Built-in touchless On/Off control
- Tune-able wave-length allows intraoral use with a range of curable composites
- 360° Light-arm rotation
- Milled billet aluminum Trolley
- Industry exclusive igus® cable management prevents cord binding and wear

UPGRADES

- ☐ 2 high-intensity Task Lights, with independent direction and focus adjustment, and wireless remote
- ☐ Adjustable VESA 100 Monitor Bracket (monitor not included)



ergonomic products
Performance Without Compromise™



www.ergonomic-products.com • 866-374-6487

Excerpted from EP Catalog:

950-075revM



Dimensions & Specifications

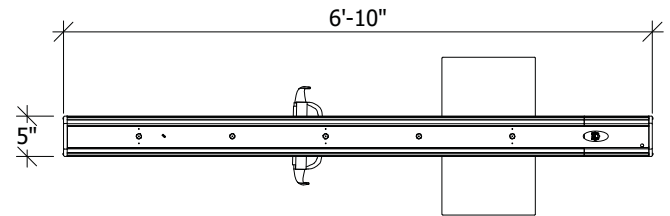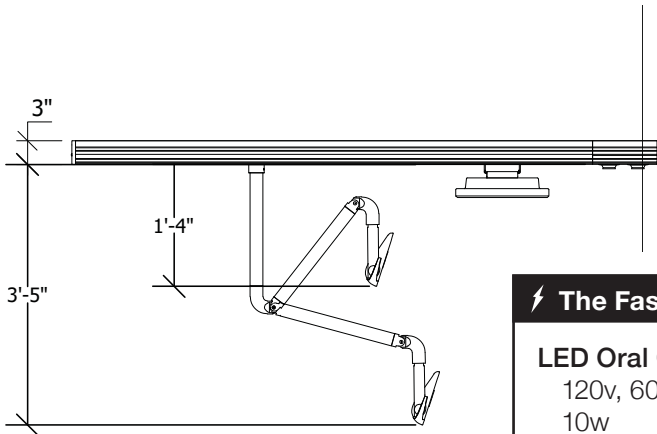

⚡ The FastTrack LS PLUS has 3 individual circuits as listed below:

LED Oral Cavity Light:

120v, 60hz
10w
Color Temp: 3500k and 5500k
Dimmable: No

LED Task Lighting (2):


120v, 60hz
8w, 3-4vdc, 680mA
Color Temp: 4000k
Dimmable: No

TV Outlet:

120v, 60hz
300w

The FastTrack PLUS is intended for mounting on a solid, flat ceiling. The Hinge-Hanger system includes a mounting-board which may be quickly affixed to blocking allowing single-person installation.

Ideal Room Location

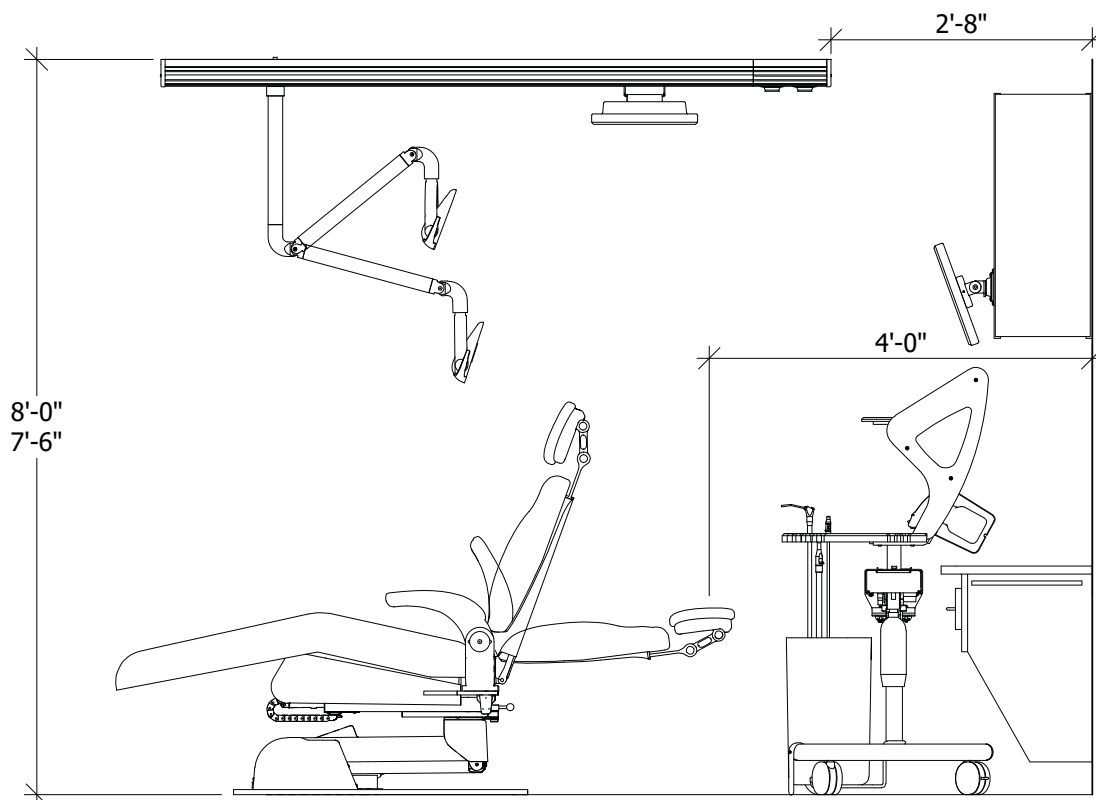

The FastTrack PLUS is intended for center-line mounting in the operator.

The layout above is represents a typical EP Workstation-equipped operator. Headwall-to-Track distances may vary in other environments. If so, position light based on reclined patient chair head-rest as indicated.