


TM

ergonomic productsTM
Performance Without Compromise[®]

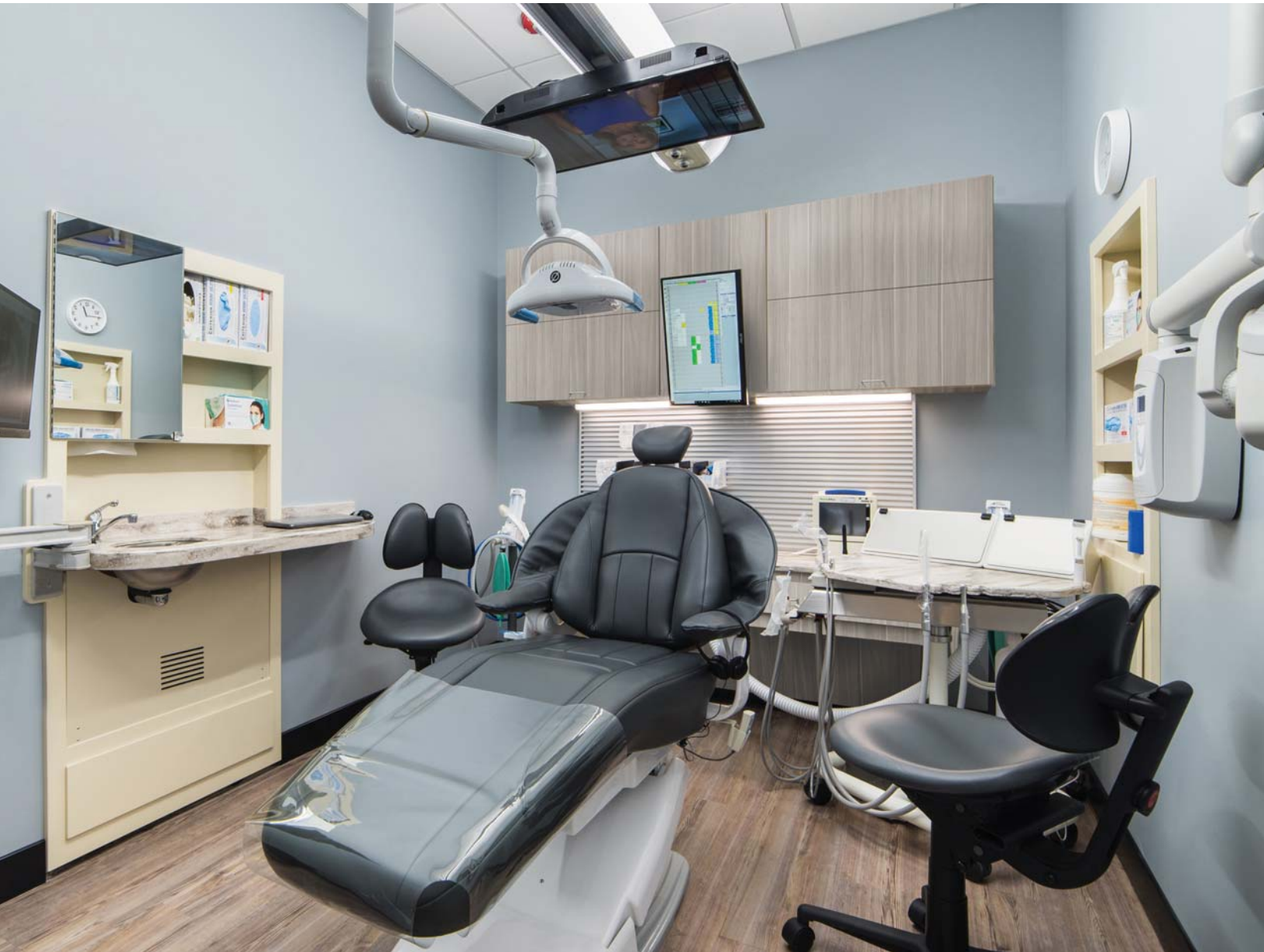

Manhattan Cabinet SystemTM

The following pages excerpted from Ergonomic Products catalog: 950-075revM

[Click here to download full catalog](#)

www.ergonomic-products.com

Manhattan Cabinet System™ (Items C1C1, M3D1 & C2C1)



Shown with options. Monitor & accessories not included.

Present patients with a **clean, welcoming** look upon entry to the operatory – without sacrificing **productivity and efficiency**.

Standard Features

C1C1 - Manhattan Upper Cabinet

- PullDown Shelf with Corian Work Surface and Keyboard Slide acts as a second tier for the Universal Workstation
- Second hydraulic PullDown
- Monitor Mounting Bracket (monitor not included)

M3D1 - Slat Wall

- Over 11 sq ft of Slat Wall to display and store accessory materials.

C2C1 - Manhattan Lower Cabinet

- Over 7 sq ft of solid Corian Work Surface
- 4 Storage Drawers
- Quick Connect Umbilical Port (standard)
- Utility Access Panel

CONFIGURATION, UPGRADE and COLOR CHOICES

PullDown Shelving Configuration and Options:

- ☐ Corian PullDown on RIGHT / Basic PullDown on left
- ☐ - or - Corian PullDown on LEFT / Basic PullDown on right
- ☐ - or - Corian PullDowns on both sides. Upgrade, Item #160-120.
(This option is needed for complete Left/Right convertibility.)

Cabinet Color/Finish Selection:

- ☐ Cabinet available in standard EP color choices.
- ☐ - or - Upgrade to woodgrain finish. (Standard EP selections)

Corian Work Surface:

- ☐ Corian available in standard and premium EP Corian color choices.

Utility Configuration

- ☐ Standard Umbilical Quick Connect
- ☐ - or - Utility Plate and Clamp-mount for 3rd party utilities

☐ Optional Sink Upgrade for C2C1 Lower Cabinet

Sink option includes stainless steel sink & chrome faucet.

The Manhattan Cabinet above is shown with a range of Slat-Wall accessories. Ergonomic Products offers the items below, and others are available from a variety of resellers.


- ☐ Med-waste Bracket
- ☐ Remote Control Holder
- ☐ Sanitation Wipes Holder
- ☐ Pad & Pen Holder
- ☐ "Goodie Bag" Brackets
- ☐ Printer Shelf
- ☐ Monitor Bracket
- ☐ Single or Triple Glove Box Holders

See page 116 for details.

TREATMENT HEADWALL



Dimensions & Utilities

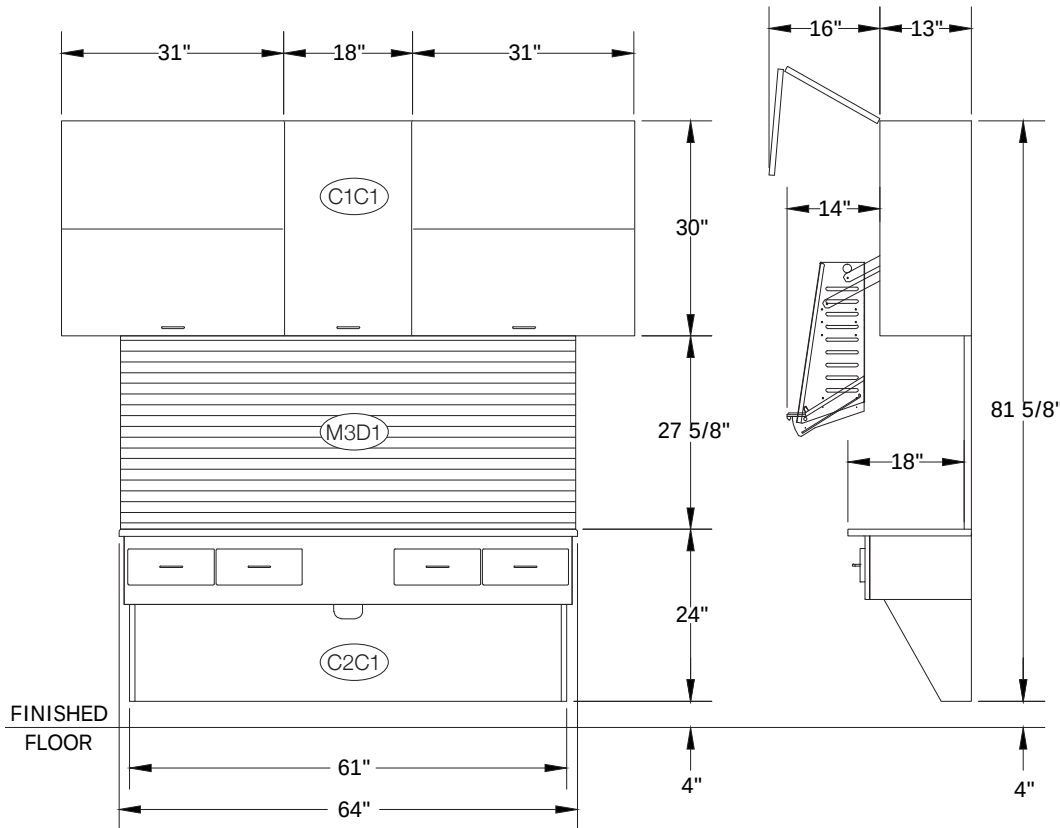

Electrical Specifications

120v AC, 60HZ, 8A
Cat 5 (or better) Data Jack

Air / Vacuum Specifications

Air: 1/4" Dental Style Female
QD, 80 psi
Vacuum: 1/2" Copper
or PVC Pipe

Plumbing Specifications

Dental Line feed:
- 1/2" Cold or Dedicated
Dental supply line with
straight 1/4 turn-stop and
3/8" compression

If C2C1 is purchased with Sink -


Sink feed:
- 1/2" Hot & Cold domestic
supply lines with straight
1/4 turn-stop and 3/8"
compression
- 1-1/2" DWV with P-Trap

Ideal Room Location

The Manhattan Cabinet System is most effective when used in a single-entry / toe-entry operatory, in conjunction with a single-tier Universal Workstation.

Used this way, much of the supplies and materials associated with treatment are removed from the patient's view on entry.

If ordered with a second Corian PullDown Shelf, the unit may be quickly converted from right-handed to left-handed use.

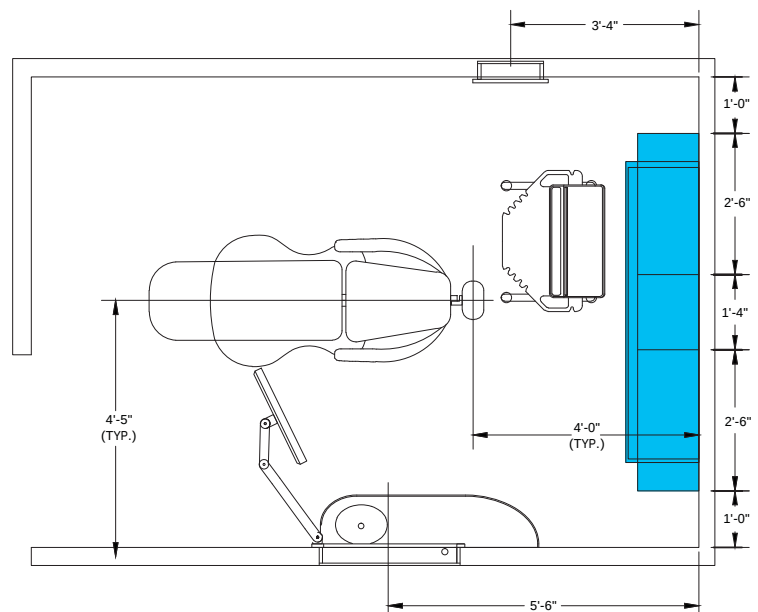

ergonomic products
Performance Without Compromise™

www.ergonomic-products.com • 866-374-6487

Excerpted from EP Catalog:

950-075revM