

# SlimLine Manhattan Cabinet System™ (Items C3B1, C3B3 & C3B4)



A **stylish look** with **enhanced efficiency** – especially for surgical operatories with space constraints.

Featuring our SlimLine aesthetic for a modern look and increased space-savings, the SlimLine Manhattan also introduces our new linear PullDown, which provides a second tier for the Workstation, without extending into the room.

## Standard Features

### C3B4 - SlimLine Manhattan PullDown

- Linear PullDown Shelf with Corian Work Surface and Keyboard Slide acts as a second tier for the Universal Workstation
- Converts from right-hand use to left-hand use in seconds

### C3B3 - SlimLine Manhattan Storage/Printer

- Center Cabinet provides storage, a Printer Shelf, and a Monitor Bracket.

### C3B1 - SlimLine Manhattan Lower Cabinet

- Nearly 5 sq ft of solid Corian Work Surface
- Open Shelf suitable for Zirc Tub placement
- Quick Connect Umbilical Port (standard)
- Utility Access Panel

Shown with options, including nitrous cutout & panel. Call for details.  
Monitor not included.

## CONFIGURATION, UPGRADE and COLOR CHOICES



### Cabinet Color/Finish Selection:

- ☐ Cabinet available in standard EP color choices.
- ☐ - or - Upgrade to woodgrain finish. (Standard EP selections)

### Corian Work Surface:

- ☐ Corian available in standard and premium EP Corian color choices.

### Utility Configuration

- ☐ Standard Umbilical Quick Connect
- ☐ - or - Utility Plate and Clamp-mount for 3rd party utilities


### Upper Cabinet Configuration

- ☐ Central Cabinet w/Printer Shelf (item C3B3 as shown)
- ☐ - or - No Printer Shelf on Central Cabinet (item C3B2)





## Dimensions & Utilities

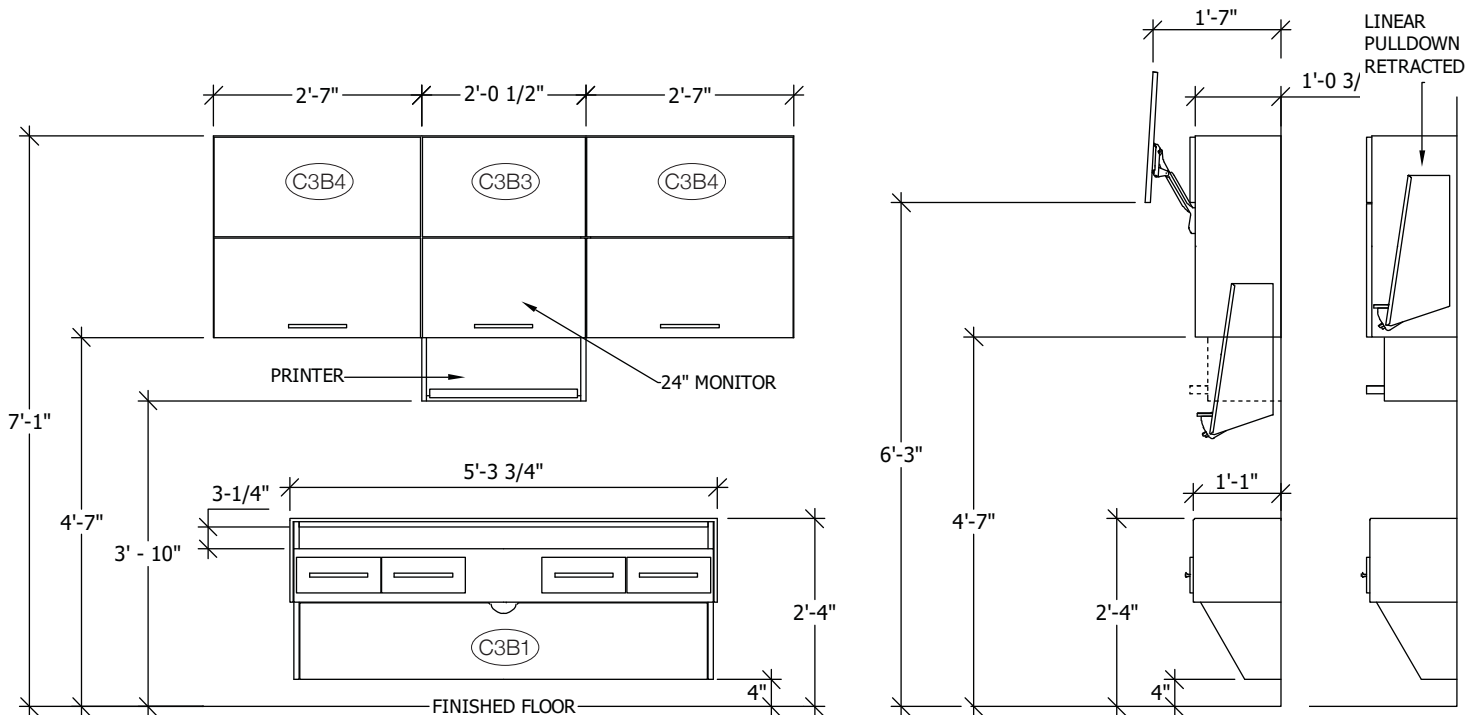

Shown with standard Printer Shelf Cabinet (C3B3). This Cabinet is also available without a Printer Shelf (C3B2).

### Electrical Specifications

120v AC, 60HZ, 8A  
Cat 5 (or better) Data Jack

### Air / Vacuum Specifications

Air: 1/4" Dental Style Female QD, 80 psi  
Vacuum: 1/2" Copper or PVC Pipe

### Plumbing Specifications

Dental Line feed: 1/2" Cold or Dedicated  
Dental supply line with straight  
1/4 turn-stop and 3/8" compression

## Ideal Room Location

Similar to our standard Manhattan Cabinet System, this SlimLine version is ideal for single-entry / toe-entry operatories and is intended to be used with a single-tier Universal Workstation.

To enhance the space-saving nature of our SlimLine products, we've modified the Corian shelf that acts as a second tier to the workstation. This variation of our Manhattan PullDown shelving travels straight up-and-down, rather than forward into the room.

With its shortened lower cabinet Work Surface, the SlimLine Manhattan can be positioned significantly closer to the head of the patient chair, reducing the overall length required for your operatories.

